

Simultaneous and accurate measurement of the dielectric constant at many frequencies spanning a wide range

Roberto Pérez-Aparicio,¹ Caroline Crauste-Thibierge,^{1, a)} Denis Cottinet,¹ Marius Tanase,¹ Pascal Metz,¹ Ludovic Bellon,¹ Antoine Naert,¹ and Sergio Ciliberto^{1, b)}
Laboratoire de Physique, Ecole Normale Supérieure de Lyon (ENS/CNRS, UMR 5672), 46 allée d'Italie, F69007 Lyon, France

(Dated: 29 January 2015)

We present an innovative technique which allows the simultaneous measurement of the dielectric constant of a material at many frequencies, spanning a four orders of magnitude range chosen between 10^{-2} Hz and 10^4 Hz. The sensitivity and accuracy are comparable to those obtained using standard single frequency techniques. The technique is based on three new and simple features: a) the precise real time correction of the amplification of a current amplifier; b) the specific shape of the excitation signal and its frequency spectrum; and c) the precise synchronization between the generation of the excitation signal and the acquisition of the dielectric response signal. This technique is useful in the case of relatively fast dynamical measurements when the knowledge of the time evolution of the dielectric constant is needed.

^{a)}Corresponding author E-mail: caroline.crauste@ens-lyon.fr

^{b)}Corresponding author E-mail: sergio.ciliberto@ens-lyon.fr

I. INTRODUCTION

Dielectric spectroscopy has provided many informations to the investigation of complex materials, and especially to the field of glass transition phenomenon (glass-forming systems) during the last decades. This technique gives direct access to the polarization of molecular dipoles allowing the study of relaxation processes. It can probe the molecular structure and dynamics of several materials like liquids, polymer composites, colloids, porous materials, or ferroelectric crystals¹. It covers a wide range of frequency, typically from 10^{-2} to 10^6 Hz, going up to 10^{10} Hz with coaxial techniques^{1,2}. Dielectric liquids of small molecules like the glass-forming glycerol³ or glass-forming polymers^{4,5} have been extensively studied.

Dielectric spectroscopy gives direct access to relaxation times⁶ like rheology, but it is a rotational relaxation one, not a translational. Nevertheless, it has been showed that on molecular liquids at equilibrium, these two relaxation times have the same temperature dependence⁷, at least above the glass transition. Recent progress in the glass transition were provided by direct measurement of a dynamical growing length scale⁸. This technique has always been improved, e.g., the important work done regarding the electrodes in order to measure charged liquids even at low frequency^{9,10}. Dielectric spectroscopy has also been coupled to scattering techniques such as X-rays¹¹ or neutrons¹².

Furthermore, recent theories on aging and out of equilibrium physics predict interesting low frequency properties of the response and fluctuations, either during the relaxation of the system towards equilibrium after a quench of the control parameter^{13,14}, or when the system is driven by external fields^{15,16}. In order to test these models it is important to measure the response of the system as a function of time after the beginning of the time dependent phenomenon¹⁷⁻¹⁹. Single frequency measurements will force to repeat the experiment many times for each selected frequency making the experiment extremely long and difficult. In those cases multi-frequency measurements are absolutely necessary.

The standard technique to perform multi-frequency response functions is to drive the system with a white noise signal. The main disadvantage of this approach is to require long averages before reaching reliable response measurements (see Section V). Indeed, the lower frequencies need long time sampling. To overcome this problem we present an innovative setup that allows simultaneous impedance measurements at several frequencies over a typical range of $10^{-2} - 10^4$ Hz with low noise and high sensibility. This has been obtained with a

careful choice of the input signal, of the current amplifier, and an appropriate correction, which allow us to measure high capacitances with very small losses ($>10^{10} \Omega$ is the typical order of magnitude). This means that our system is particularly useful for the study of relaxation processes in polymer films which have very high resistivity. The length of the measurement time window is imposed only by the smallest frequency that one needs to resolve because all other frequencies are measured simultaneously.

The paper is organized as follows. In Section II we give some general, basic background on the dielectric spectroscopy. Section III describes the experimental principle, IV the experimental setup, and V the measurement procedure. Application of the technique to the study of some polymer films is presented in Section VI. We conclude in Section VII.

II. BACKGROUND

Dielectric spectroscopy investigates the frequency dependent dielectric properties of a material¹. It is based on the interaction of an external field with the electric dipole moment of the sample, often expressed by the complex dielectric permittivity or dielectric constant:

$$\varepsilon(\omega) = \frac{D(\omega)}{E(\omega)}, \quad (1)$$

where $D(\omega)$ is the dielectric displacement and $E(\omega)$ the electric field, both dependent on the angular frequency ω . $\varepsilon(\omega)$ can be expressed as

$$\varepsilon(\omega) = \varepsilon'(\omega) - i\varepsilon''(\omega), \quad (2)$$

being $\varepsilon'(\omega)$ the stored permittivity, which is related to the stored energy by the sample, and $\varepsilon''(\omega)$ the loss permittivity, which indicates the loss energy. Their loss tangent is

$$\tan(\delta) = \frac{\varepsilon''(\omega)}{\varepsilon'(\omega)} \quad (3)$$

where δ is the phase shift between $D(\omega)$ and $E(\omega)$. The dielectric constant ε depends on the dielectric material, the frequency, the temperature, age.

The orientation polarization of dipoles due to applied disturbances allows the investigation of the dielectric relaxation of the material in the frequency range below 10^9 Hz. Moreover, if those dipoles are of molecular origin, the dielectric spectroscopy allows to investigate molecular motions. The experimental techniques and procedures for dielectric spectroscopy

depend on the frequency range of interest. At low and intermediate frequencies ($10^{-3} - 10^7$ Hz) the sample can be modeled as a parallel plate capacitor where the sample is the dielectric insulator in between the two electrodes²⁰. Assuming no edge capacitances, the capacitance of an ideal parallel plate capacitor reads

$$C(\omega) = \frac{\varepsilon_0 \varepsilon'(\omega) S}{d}, \quad (4)$$

where d is the distance between electrodes, S the surface occupied by the sample and ε_0 the vacuum permittivity. Using Eq. 2 in Eq. 4, we see that the impedance $Z(\omega)$ of the sample is modeled by a resistor $R(\omega)$ in parallel with a capacitor $C(\omega)$ and it can be expressed as :

$$\frac{1}{Z(\omega)} = \frac{1}{R(\omega)} + iC(\omega)\omega \quad (5)$$

with $\varepsilon'(\omega) = C(\omega)(\varepsilon_0 S)^{-1}d$ and $\varepsilon''(\omega) = (\varepsilon_0 S R(\omega)\omega)^{-1}d$.

The complex impedance $Z(\omega)$ can be measured by imposing the voltage $V(\omega)e^{i\omega t}$ and measuring the current $I(\omega)e^{i\omega t}$ flowing through:

$$Z(\omega) = \frac{V(\omega)}{I(\omega)} \quad (6)$$

As ε'' is low on pure dielectric systems this means that the resistive part of the impedance is quite high and the in-phase current becomes hardly detectable.

III. THE PRINCIPLE OF THE MEASURE

So far we described the ideal capacitor, but in the real experiments there are other parameters to take into account. On the one hand, corrections due to the measurement system like the influence of the extra capacitance in the cables. On the other hand, sample preparation must be very careful to avoid difficulties, such as sample geometry, bad electric contacts between electrodes and sample, edge capacitances, or electrode polarization. Sample geometry has to be well known. We do not discuss these standard experimental details but we focus on the description of the features of the electronics, the driving signal, and the data analysis that we use in our experimental setup.

Our goal is to measure $\varepsilon(\omega)$ simultaneously in a frequency range covering 4 orders of magnitude in ω . As dielectric losses are in general small in pure dielectric (such as polymers), any spurious effect must be taken into account in order to have a good accuracy on the wide

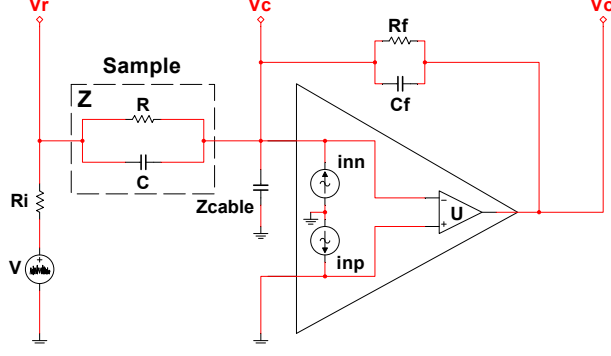


FIG. 1. Scheme of the working principle of the main amplifier of our device.

frequency range. To define notation and to clarify the principle of our technique let us recall basic electronic concepts. In order to measure $I(\omega)$ we use a current amplifier whose entrance impedance is just the sample impedance $Z(\omega)$ as sketched in Fig. 1. In order to analyze the drawbacks of this scheme let us assume first that the amplifier U in Fig. 1 is an ideal one. Thus, we suppose that the bias currents at both entrances are zero ($i_{nn} = i_{np} = 0$), and that the open-loop amplification of U is infinite for all frequencies. This implies that $V_c = 0$ because the positive entrance is at the ground potential. Therefore, the impedance $Z_{cable} \rightarrow \infty$ (given by the parallel of the capacitances of the cables and of the amplifier entrance) does not play any role in this case. From the standard circuit analysis, the current balance between the current flowing in the feedback impedance Z_f ($C_f // R_f$) and the entrance current we obtain

$$\frac{V_o}{Z_f} = -\frac{V_r}{Z}, \quad (7)$$

where V_r is the reference and V_o the output potentials (see Fig. 1). We have omitted the obvious dependence on ω . Thus, the impedance Z can be calculated as

$$Z = -Z_f \frac{V_r}{V_o} \quad (8)$$

Nevertheless, considering a real amplifier, (i.e. with a finite open loop gain) one has to take into account in the current balance (Fig. 1) that $V_c \neq 0$. Thus, Z_{cable} cannot be neglected and the current balance is :

$$\frac{V_o - V_c}{Z_f} + \frac{V_r - V_c}{Z} - \frac{V_c}{Z_{cable}} + \delta I_o = 0, \quad (9)$$

where we have considered that the random output current δI_o takes generically into account all the sources of noise produced by the entrance bias currents, the voltage noise, and the

component thermal noise. In order to get Z we first multiply Eq. 9 by the conjugate complex $V_r^*(\omega)$ of $V_r(\omega)$ and we take an average over several realizations:

$$\frac{\langle V_o V_r^* \rangle - \langle V_c V_r^* \rangle}{Z_f} + \frac{|V_r|^2 - \langle V_c V_r^* \rangle}{Z} - \frac{\langle V_c V_r^* \rangle}{Z_{cable}} = 0, \quad (10)$$

where we took into account that δI_o is uncorrelated with V_r thus $\langle \delta I_o V_r^* \rangle = 0$. We consider the transfer functions of both the signal and the correction:

$$G_{sig}(\omega) = \frac{\langle V_o(\omega) V_r^*(\omega) \rangle}{|V_r|^2(\omega)} \quad (11)$$

$$G_{corr}(\omega) = \frac{\langle V_c(\omega) V_r^*(\omega) \rangle}{|V_r|^2(\omega)} \quad (12)$$

Using Eq. 10, and the definitions of G_{sig} and G_{corr} , the impedance of the sample can be calculated as

$$\frac{1}{Z} = \left(\frac{G_{corr}}{Z_{cable}} + \frac{G_{corr} - G_{sig}}{Z_f} \right) \frac{1}{1 - G_{corr}} \quad (13)$$

The knowledge of V_c is very important in order to have reliable estimations of the dielectric constant at relatively high frequencies when the amplifier open loop gain begin to decrease and the phase shifts are not negligible with respect to the one that we want to measure. In our experimental setup we actually measure V_c as we describe in the next section.

IV. EXPERIMENTAL SETUP

The real scheme of our device is presented in Fig. 2. The characteristics of each component are presented in Table I. The source signal $V(t)$ is sent to a first operational amplifier U_1 . This amplifier is a voltage follower or a non-inverting buffer. This leads to decouple the impedance of the source and the sample with a negligible additional noise. For that, we choose the LME 49990 amplifier, which has very low voltage noise. The output of U_1 is sent to the sample via R_2 which is the equivalent to R_i in Fig. 1. Note that $R_i = R_2$ is very useful in order to avoid self oscillations of the operational amplifier U_1 when large C are used (see Appendix). The signal sent to the sample is V_r as defined in the previous section.

V_r is also sent to a secondary buffer amplifier U_2 (LME 49990), and the output of this amplifier is used as the reference signal V_{ref} (assuming $V_{ref} = V_r$ because the phase-shift of U_2 is negligible in our frequency range).

TABLE I. Technical details of the components of the circuit of Fig. 2. U are the different amplifiers of types LME 49990 (low voltage noise, $1.4 \text{ nV}/\sqrt{\text{Hz}}$ at 10 Hz) or AD549 (low current noise, $0.5 \cdot 10^{-15} \text{ A}_{rms}$, $f = 0.1$ to 10 Hz, bias current 30 fA); R are resistances.

U	R (Ω)	
U ₁ LME 49990	R ₁	10^4
U ₂ LME 49990	R ₂ = R ₅ = R ₇ = R ₁₀ = R ₁₄ = R _i	10^3
U ₃ AD549	R ₃ = R ₉ = R ₁₁	50
U ₄ LME 49990	R ₄ = R _f	10^6
U ₅ AD549	R ₆ = R ₈ = R ₁₃ = R ₁₅	$4 \cdot 10^3$

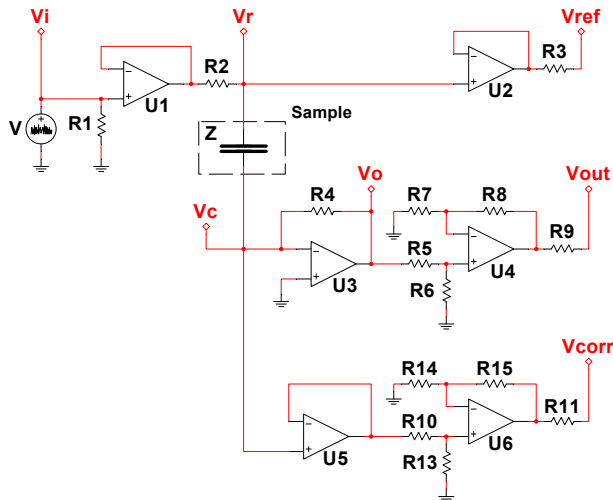


FIG. 2. Scheme of our real circuit. U₁, U₂, and U₄ are low noise voltage amplifiers. U₃ and U₅ are low noise current amplifiers. Details of the components are on Table I.

The current flowing through the sample can be very small if the sample impedance is very high so this part of the circuit uses amplifiers with a low current noise and bias current (AD549). There are two acquisition paths: in the main one, the voltage from the sample goes through a first amplifier U₃ with $R_4 = R_f = 10^6 \Omega$. This amplifier is the main one of the circuit and is equivalent to the one modeled in Fig. 1 and by Eq. 9. It is followed by U₄ (LME 49990) which amplifies the voltage by a factor 4 taken into account in data treatment. The second path is the measurement of the negative entrance potential of amplifier U₃, that is to say V_c in Eq. 9. This signal also goes through a buffer amplifier U₅ and is amplified 4

times by U_6 , getting V_{corr} (i.e. $V_{corr}/4 = V_c$).

In order to accurately measure impedances, we have to take into account all the contributions from the different components in the electronic device. The impedance of the cables Z_{cable} between the sample and the amplifiers has to be considered

$$Z_{cable} = \frac{1}{i\omega(C_{cable} + C_e)}, \quad (14)$$

where $C_{cable} = 50 \cdot 10^{-12}$ F is the capacitance of 0.5 m long cables we used. The entrance capacitance $C_e \simeq 10^{-12}$ F of U_3 can be neglected. We made careful choice of good quality cables between the sample and the amplifier –other cables have less sensitive contribution to the noise in the circuit– in order to reduce noise in the signal. We also have to take into account the feedback impedance of the amplifiers

$$Z_f = \frac{1}{R_f^{-1} + i\omega C_f}, \quad (15)$$

where $R_f = 10^6 \Omega$ is the feedback resistance and $C_f = 10^{-13}$ F its parallel capacitance (Fig. 1). This capacitance is negligible compared to the impedance of R_f in our frequency range.

The signals V_{ref} , V_{corr} , and V_{out} are sent to a National Instrument PXI-4472 acquisition card (24 bits resolution). Data acquisition is performed via a Labview home-made program and data treatment using Matlab to compute the transfer functions $G_{sig}(\omega)$ (Eq. 11) and $G_{corr}(\omega)$ (Eq. 12) using standard FFT algorithm. The transfer functions of the three amplifiers composed by U_2 , U_4 , and (U_2+U_6) are measured in order to correct the phase of the three measured signals.

V. THE MEASUREMENT PROCEDURE

In order to perform a simultaneous and precise measurement of Z in a wide frequency range using the device described in the previous section we have optimized the frequency dependence of the signal V_r . Indeed the standard technique employed to perform a response measurement on a wide band is to excite the system with an almost white noise in the band of interest and to measure simultaneously the input and output signals. Then, by taking FFT of these two signals input and output, one can compute the transfer function G_{sig} by a very long average of the response function on many configurations of the noise. This system

although very practical and adaptable on almost all cases presents several drawbacks in the case of a circuit like the one presented in the previous sections. The first is that in the circuit of Fig. 2 the high frequencies are amplified much more than the low ones. Therefore, one has to compensate for this effect. The second drawback is using a random noise at the input, its phase is rapidly changing and many averages are required in order to smooth this phase noise and to have a reliable measurement. In order to avoid these difficulties and reduce the number of averages in our system we apply the following method.

A. Method

We first fix a measurement interval T and a sampling frequency $f_s = N/T$, where N is an integer. Then, we generate a signal composed by a sum of sinusoids whose pulsations ω_k are harmonics of the slowest frequency $2\pi/T$, specifically

$$V(t) = \sum_k V_{0,k} \sin \omega_k t, \quad (16)$$

where $V_{0,k}$ are the amplitudes in volts. The amplitudes are chosen to be $V_{0,k} \propto 1/\omega_k$ because favoring low frequencies allows us to compensate the increasing gain of the amplifier at high frequencies. These kinds of signal are numerically generated and recorded in the computer memory which is then sampled by the same clock at the sampling frequency f_s , which is also used for sampling the acquired signals V_r , V_o , and V_c . In such a way all phase noises are suppressed as each frequency contains a fix integer number of acquired samples. We use two different devices to generate the signal either a National Instrument PXI 5411-card generator or an Agilent 33500B Series waveform generator. As an example we plot in Fig. 3 the spectra of V_{ref} , V_{out} , and V_{corr} signals obtained using $C = 10^{-8}$ F. In this specific example we select twenty-five frequencies $f_k = \omega_k/(2\pi)$ in a logarithmic scale between a minimum value of 0.125 Hz and 10^3 Hz. The sampling rate is $2^{15} = 32768$ sample/s. The peak of the noise at about 3 kHz is explained in the Appendix. Notice that the amplitude of the correction signal is not negligible above 100 Hz (Fig. 3c) and it is important to keep it into account in the calculation of Z . To summarize, this electronics combined with the described data analysis can measure with a good accuracy $R \leq 10^{12}$ Ω and $C > 10^{-11}$ F.

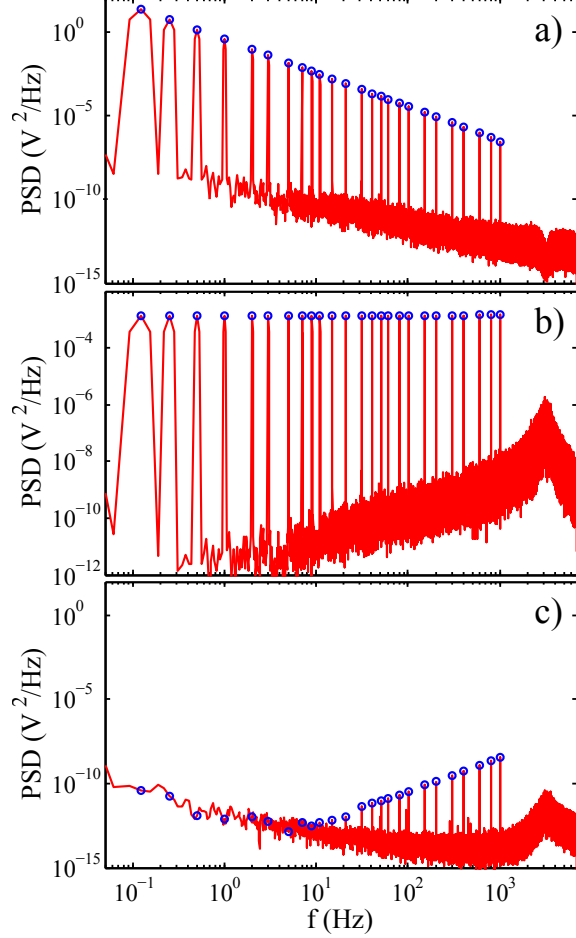


FIG. 3. Power spectral density of (a) V_{ref} , (b) V_{out} , and (c) V_{corr} signals. Circles point to the chosen frequencies.

B. Data Analysis

From the measured impedance (Eq. 13) we can calculate the resistance and the capacitance of the dielectric as

$$R = \left[\text{Re} \left(\frac{1}{Z} \right) \right]^{-1} \quad (17)$$

$$C = \text{Im} \left(\frac{1}{Z} \right) \omega^{-1} \quad (18)$$

In order to check the efficiency of our device in the chosen frequency window, we use several impedances Z whose values have been checked by standard apparatus²¹. In Fig.4a we plot the measured values of the capacitances and in Fig.4 b) those of the measured resistances. This plot illustrates the whole dynamics of our system ranging from 10^{-12} F (at 10^{-13} F the

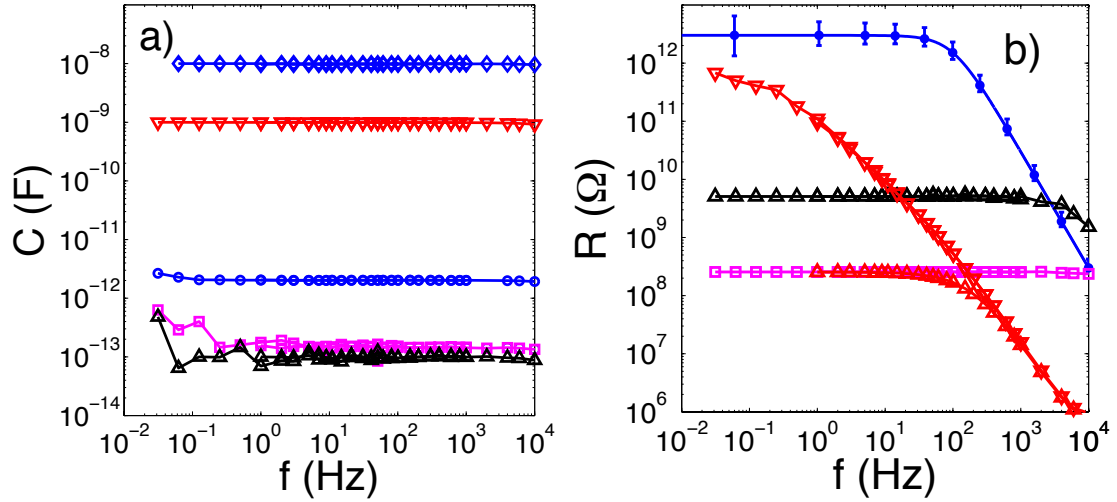


FIG. 4. Capacitances (a) and resistances (b) of several components used to check the accuracy and the limits of our apparatus. Three capacitances of 2.5 pF (●), 1 nF (▽) and 10 nF (◇) are plotted. The losses of the 2.5 pF (● with error bars) and 1nF (▽) capacitances are plotted in (b), whereas those of the 10 nF in fig.5b. The capacitance of 2.5 pF has very small losses and this allows us to fix the limits and the losses of our apparatus (blue line with error bars in b). In (b) we also plot the value of two pure resistances of 250 MΩ (□) and 5 GΩ (△). These two resistances have of course a small parallel capacitance of about 0.1 pF shown in (a) (□) and (△), which also fixes the accuracy of our system. In (b) we also plot the measured resistance (△) when Z is the parallel of the 1 nF capacitance with the 250 MΩ resistance. Notice that for each Z the whole frequency interval is covered by only two simultaneous measurements in two four order of magnitude intervals: $[10^{-2} - 10^2 \text{ Hz}]$ and $[1 - 10^4 \text{ Hz}]$

error is big but it can be still used) to 10 nF and larger if a smaller R_f is used²². The 2.5 pF capacitance is an air capacitance with extremely small losses and this allows us to determine the losses of our circuit as a function of frequency which are plotted in Fig.4b. In order to appreciate more precisely the accuracy we plot in Fig.5a) the measured values of the 10 nF capacitance and Fig.5b) its losses in the frequency range 0.1 Hz – 1 kHz. We see that the fluctuations of the measurement are less than 1/1000 in all the frequency range. In Fig.5b we compare the losses of the 10 nF capacitance with a measured resistance $R_e = 5 \text{ G}\Omega$. The capacitance C has an almost constant value, whereas its losses increase, which correspond to

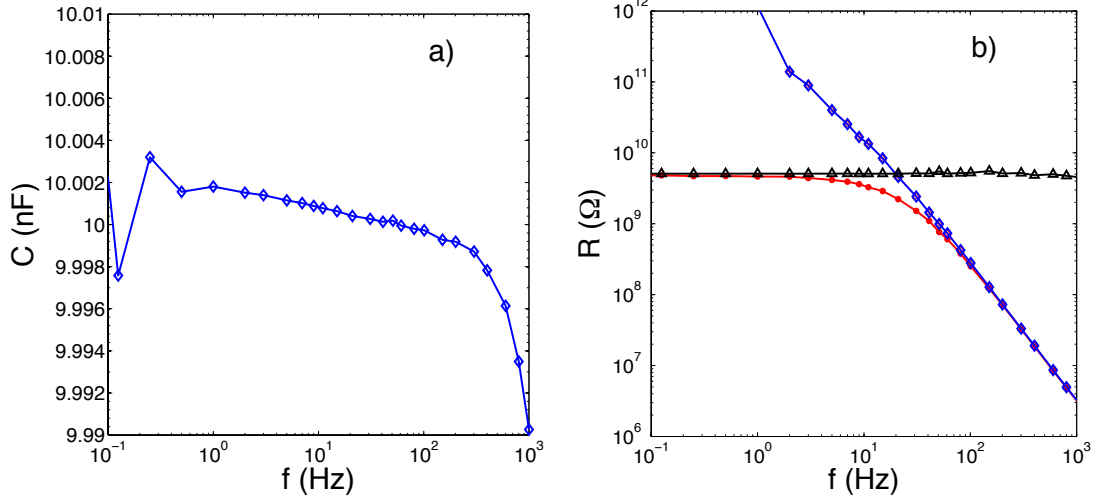


FIG. 5. (a) Expanded view of the measured 10 nF capacitance of Fig.4a. Its losses are plotted in (b) together with the pure resistance of 5 G Ω (Δ). We also plot the measured total resistance \bullet of Z composed by the parallel of the 10 nF with the 5 G Ω resistance. The dashed line represent the results of the measurement without using the correction signal V_c in the calculation. The whole measurement in the range 0.1 Hz-1000 Hz takes less than 60 s (i.e. 8 periods of the smallest frequency).

a decrease of $R(\omega)$ as a function of ω (see eq.5 and Fig.5b). When C is mounted in parallel to R_e then the capacitance C is unchanged whereas the measured resistance is now the parallel of $R_e//R(\omega)$ which is constant at low frequencies up to around 10 Hz, and then, decreases as a function of the increasing frequency due to the capacitive losses²¹. It is important to stress that, using our method, we need to average less than 60 s, in order to measure simultaneously the values of R and C on 4 orders of magnitude in frequency, plotted in Fig. 5. As already mentioned, standard techniques use a white noise in the bandwidth of interest. By applying such a signal to our circuit, the result is much more noisy (even impossible to find correct values) than that obtained using our specifically chosen input signal described in the previous section (see Eq.16). In Fig. 5 we can also see the importance of the correction signal V_c (Fig. 3c in order to properly calculate the values of R and C in the whole chosen bandwidth. Finally we stress that because of the limited memory of our voltage generator we always measure simultaneously at most 4 order of magnitude in frequency in the range 0.01 Hz – 10 kHz. Furthermore the accuracy is not constant as a function of frequency. It

is 5 % below 0.1 Hz and above 1 kHz about 1 % in the range 0.1 Hz – 1 kHz

VI. MATERIALS AND APPLICATIONS

In the case of a pure capacitance, from the values of R and C we can calculate both ε' and ε'' :

$$\varepsilon' = \frac{C}{C_0} \quad (19)$$

$$\varepsilon'' = \frac{1}{R\omega C_0} \quad (20)$$

with $C_0 = S\varepsilon_0 d^{-1}$. Therefore, $\tan(\delta)$ (Eq. 3) can be calculated as

$$\tan(\delta) = \frac{1}{RC\omega}, \quad (21)$$

which is independent of the geometry of the sample.

We check our experimental device with well-known polymer films. As already mentioned, our experimental device allows to perform simultaneous measurement of impedance at several frequencies in a range from 10^{-2} to 10^3 Hz with low noise and high sensitivity. First of all, we investigate an extrusion film MAKROFOL[®] DE 1-1 000000 (from BAYER) based on Makrolon[®] polycarbonate (PC) of 125 μm of thickness and with a glass transition temperature (T_g) of about 150 °C. The sample consist of a PC film of 18 cm of diameter. We measure the impedance and compute C and R from 22 to 167 °C. Fig. 6 shows R and C at 22 °C. R presents huge values at low frequencies decreasing with frequency. C presents a small dependence in frequency. From C we calculate the dielectric constant ε (Eq. 4) to be similar to the value found in the literature²³ ($\varepsilon(60 \text{ Hz}) = 3.17$). $\tan \delta$ (Eq. 21) as a function of temperature is presented in Fig. 7. A peak corresponding to the segmental relaxation of PC (glass transition) appears from low to high frequencies when increasing temperature (well-known in polymers, see e.g. ref. 2). The advantage here is that one can do easily multi-frequency dynamical measurements changing the rate of the temperature ramp.

Then we investigate a film of poly(vinyl-acetate) (PVAc) in order to study the aging phenomenon. The PVAc ($T_g \sim 34^\circ\text{C}$) is cycled between 19 and 53°C, with a very fast thermal quench (around 30°C/min). The electrodes are made out of thin films of aluminum (150 nm) deposited by evaporation over a sapphire disc. Such crystal plate is chosen for its

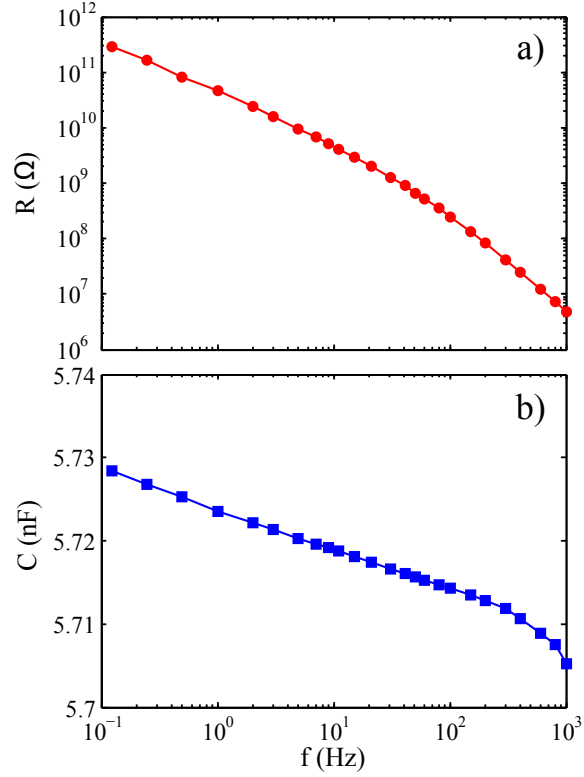


FIG. 6. Resistance (a) and capacitance (b) of a PC film of $125 \mu\text{m}$ of thickness and 18 cm of diameter as a function of frequency at room temperature.

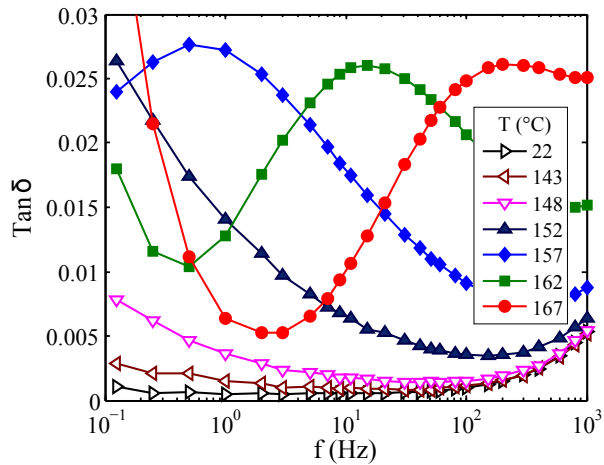


FIG. 7. $\text{Tan } \delta$ of a PC film of $125 \mu\text{m}$ of thickness and 18 cm of diameter as a function of frequency at different temperatures.

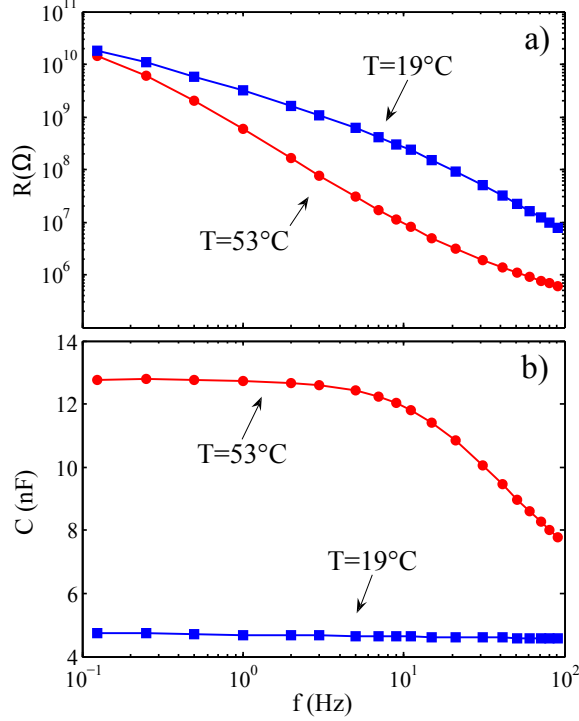


FIG. 8. Resistance (a) and capacitance (b) of a PVAc sample as a function of frequency at two different temperatures.

high thermal conductivity because heating and cooling are done through the substrate. The polymer is dissolved in chloroform, then poured over the lower electrode. After careful drying of the solvent, the second aluminum thin film is deposited over the dielectric. Connexions are fixed on the electrodes thanks to conductive glue. The thickness of the dielectric is around $1 \mu\text{m}$, such that the overall capacitance is a few nF in order to match the amplifier's impedance. Fig. 8 shows resistance and capacitance of PVAc above and below T_g as a function of frequency. The α -relaxation appears in the measurement range around 53°C but not at lower temperature (19°C). Fig. 9 shows the time evolution of R and C after a fast thermal quench from 53°C to 19°C at 1, 11, and 91 Hz. The relaxation of R and C are due to the aging of the polymer below the glass transition going toward equilibrium. This technique allows to work on a time range limited at small times by the speed of the quench (typically one minute). Our device is absolutely necessary in this kind of dynamical studies.

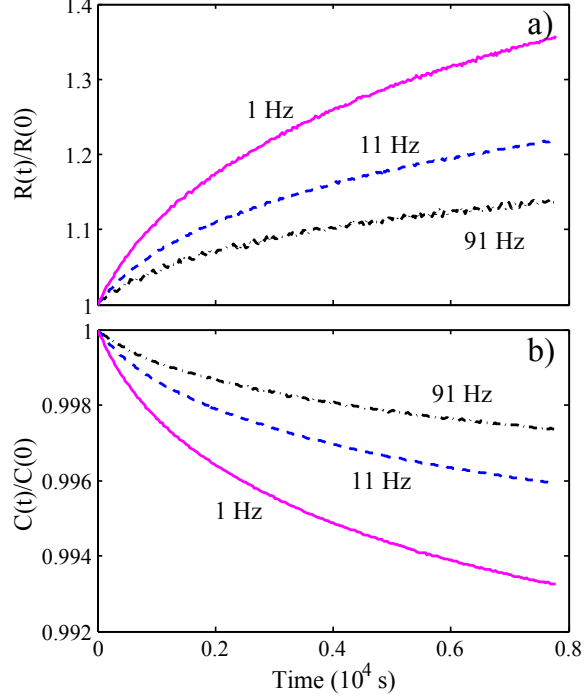


FIG. 9. Resistance (a) and capacitance (b) of a PVAc sample as a function of the aging time after a quench from 53 to 19 °C at 1, 11, and 91 Hz. Both magnitudes are normalized to their values just after the quench (time = 0): $1.95 \cdot 10^9 \Omega$ (1 Hz), $2.5 \cdot 10^8 \Omega$ (11 Hz), and $3.4 \cdot 10^7 \Omega$ (91 Hz) for the resistance; $5.99 \cdot 10^{-9} \text{ F}$ (1 Hz), $5.89 \cdot 10^{-9} \text{ F}$ (11 Hz), and $5.82 \cdot 10^{-9} \text{ F}$ (91 Hz) for the capacitance.

VII. CONCLUSION

We present here a technique for impedance measurements, specially designed for highly resistive samples and time dependent phenomena. Indeed it allows us to measure the complex dielectric constant simultaneously in a wide frequency range. The electronics combined with the described data analysis allow to measure with a good accuracy $R \leq 10^{12} \Omega$ and $C > 10^{-11} \text{ F}$.

This technique is based on two main elements. First of all, a home-made amplifier specially built to measure small phase shifts in the whole frequency range of interest. The proposed method based on a correction signal is very efficient for obtaining accurate values of both on the real and the imaginary part of the impedance, even in the cases of very small phase shift.

Secondly, the reference signal is carefully chosen to give optimal accuracy over the whole frequency range, which is typically $10^{-2} - 10^4$ Hz. We use a sum of sinus with decreasing intensity at high frequency in order to keep a constant signal-to-noise ratio with frequency. The chosen frequencies of the input signal are exact subharmonic of the sampling frequency, which is used to generate the input signal and to sample the outputs. Compared to the use of a white noise, this signal has well-defined phases, ensuring the good phase accuracy we need. Furthermore, measuring known frequencies allows shorter acquisition times than using white noise in order to get high precision on the impedances. Then, this device is particularly useful for time evolving impedances studies, especially in the field of aging of glass-forming polymers. Four orders of magnitude in frequency, chosen in the range of $10^{-2} - 10^4$ Hz, can be simultaneously investigated over typical time range from one minute to arbitrary long time. We also present here few results on two polymer samples: a PC film evolving to the glass transition with temperature, and a PVAc sample evolving with time after a fast quench. This last result demonstrates the rich field of non-equilibrium physics that can be probed using this device.

ACKNOWLEDGMENTS

R. Pérez-Aparicio acknowledges the funding by Solvay. We also thank Debjani Bagchi who made use of the setup presented here.

Appendix: Response to noise

In this appendix we give an explanation of the noise peak in the output spectrum of Fig. 3. The amplification of the AD549 is very well approximated at low frequencies by

$$A = \frac{A_o}{1 + i \omega / \omega_c} \quad (\text{A.1})$$

with $A_o \simeq 10^6$ and $\omega_c \simeq 2 \pi$ rad/s.

The analysis of the circuit of Fig. 2 shows that the spectrum of the output signal $|V_o(\omega)|^2$ is related to that at the input $|V_i(\omega)|^2$:

$$|V_o(\omega)|^2 = \frac{(R_f \omega C)^2 A_o^2 |V_i(\omega)|^2}{D(\omega)} \quad (\text{A.2})$$

where

$$D(\omega) = (1 - \omega^2(R_i + R_f)C/\omega_c + A_o)^2 + [(R_i + R_f) C + 1/\omega_c + A_o R_i C]^2 \omega^2 \quad (\text{A.3})$$

Here we have used the approximation $R \rightarrow \infty$ and $C_f \simeq 0$ (see Fig. 1). These are very good for the purpose of understanding the maximum of the noise. This spectrum has a resonance at $\omega \simeq \sqrt{A_o \omega_c / [(R_i + R_f)C]}$ which is rather broad because of the strong damping term $A_o R_i C$. Inserting the numerical value we see that the resonance is at about 3.9 kHz and the broadening a few Hz. These values are in agreement with the enhancement of the noise in the power spectrum of Fig. 3b, if we consider that the noise of the input signal, Fig. 3a, has a wide bandwidth with a slow decrease as function of frequency. In Fig. 10 we plot an output spectrum computed from Eq. A.2 using a realistic input spectrum, which is plotted in the same figure. We see that the output spectrum is very close to the one measured with $C = 10^{-8}$ F.

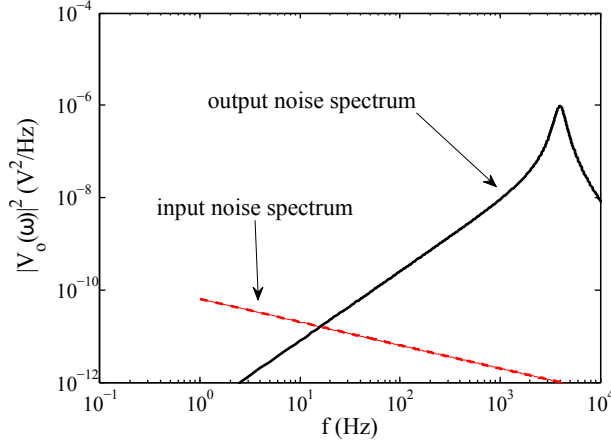


FIG. 10. Output spectrum (line) computed using in Eq. A.2 the plotted input spectrum (dashed-line).

REFERENCES

- ¹F. Kremer and A. Schonhals, *Broadband Dielectric Spectroscopy* (Springer-Verlag, Berlin, 2003).
- ²P. Lunkenheimer, U. Schneider, R. Brand, and A. Loidl, *Contemporary Physics* **41**, 15 (2000).

- ³P. Lunkenheimer, A. Pimenov, M. Dressel, Y. G. Goncharov, R. Böhmer, and A. Loidl, *Phys. Rev. Lett.* **77**, 318 (1996).
- ⁴A. J. Bur, S. C. Roth, and M. McBrearty, *Review of Scientific Instruments* **73**, 2097 (2002).
- ⁵D. Gómez, A. Alegría, A. Arbe, and J. Colmenero, *Macromolecules* **34**, 503 (2001).
- ⁶S. Havriliak and S. Negami, *Polymer* **8**, 161 (1967).
- ⁷I. Chang and H. Sillescu, *The Journal of Physical Chemistry B* **101**, 8794 (1997).
- ⁸C. Crauste-Thibierge, C. Brun, F. Ladieu, D. L'Hôte, G. Biroli, and J.-P. Bouchaud, *Phys. Rev. Lett.* **104**, 165703 (2010).
- ⁹B. L. Mellor, E. C. Cortés, S. Khadka, and B. A. Mazzeo, *Review of Scientific Instruments* **83**, 015110 (2012).
- ¹⁰P. B. Ishai, Z. Sobol, J. D. Nickels, A. L. Agapov, and A. P. Sokolov, *Review of Scientific Instruments* **83**, 083118 (2012).
- ¹¹I. Šics, A. Nogales, T. A. Ezquerra, Z. Denchev, F. J. Baltá-Calleja, A. Meyer, and R. Döhrmann, *Review of Scientific Instruments* **71**, 1733 (2000).
- ¹²M. Jiménez-Ruiz, A. Sanz, A. and Nogales, and T. A. Ezquerra, *Review of Scientific Instruments* **76**, 043901 (2005).
- ¹³L. F. Cugliandolo, J. Kurchan, and L. Peliti, *Phys. Rev. E* **55**, 3898 (1997).
- ¹⁴J.-L. Barrat and W. Kob, *Europhys. Lett.* **46**, 637 (1999).
- ¹⁵J.-L. Barrat and L. Berthier, *Phys. Rev. E* **63**, 012503 (2000).
- ¹⁶J. R. Gomez-Solano, A. Petrosyan, S. Ciliberto, R. Chetrite, and K. Gawedzki, *Phys. Rev. Lett.* **103**, 040601 (2009).
- ¹⁷L. Bellon, S. Ciliberto, and C. Laroche, *Europhys. Lett.* **53**, 511 (2001).
- ¹⁸D. Hérisson and M. Ocio, *Phys. Rev. Lett.* **88**, 257202 (2002).
- ¹⁹L. Buisson, S. Ciliberto, and A. Garcimartín, *Europhys. Lett.* **63**, 603 (2003).
- ²⁰We do not discuss here this obvious problem because it is not the subject of the paper and it can be easily taken into account in our experimental setup by a guard ring.
- ²¹The values of this test circuit have been checked both by a Novocontrol equipment and a Stanford lock-in amplifier (SR830).
- ²²This is done using $R_f = 1 \text{ M}\Omega$ but it can be extended to larger capacitance by reducing its value.
- ²³J. E. Mark, *Polymer Data Handbook* (Oxford University Press, Oxford, 1999).